



Using the Automated External Defibrillator (AED)

A First Aid Device That Guides the User Step by Step

"A device designed to help people in cases of sudden cardiac arrest, operating with audio or visual guidance until specialists arrive"



Why Do We Care About the AED?

Time is an important factor in sudden cardiac arrest cases



A Quick Response
in Critical Cases



Voice Instructions
That Guide the User
Precisely



Available in
Public Facilities



It Supports CPR
Within the Initial
Response



Community Readiness
to Handle Crises



Helping the Injured Person
Until the Specialists
Arrive

Remember: Having an AED alone isn't enough — most important is knowing its location, requesting help, and following its instructions precisely



What Is an AED?

An automated electronic device that analyzes the heart rhythm and guides the user



Audio or visual instructions



It helps step by step



It analyzes the heart's condition automatically



No shock except on the device's recommendation



Used with CPR when needed



Designed to be clear and easy

Note: This presentation is for awareness and does not replace approved hands-on training in CPR and AED



When Do We Use an AED?

Signs that call for immediate action



The person doesn't respond when called

1

Not breathing normally

2

Abnormal breathing resembling gasping

3

A sudden collapse and loss of responsiveness

4

No clear movement or response

5

Securing the area before approaching

6

Ask for an AED while calling 999

7

Urgent Action: Make sure the area is safe, call 999, and ask for an AED to be brought if available



The Basic Rule Before Using an AED

Check for hazards before approaching



Fire or smoke



**Water around
the injured
person or the
device**




**Exposed
electrical
wiring**



**Nearby vehicle
traffic**



**A risk of falling
or crowding**

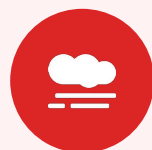

**A gas or fuel
leak**

Important Message: Don't use the device in an unsafe location, and don't put yourself or others in danger



Initial Response Steps

Act in clear order and calmly

Start CPR if needed

5

Turn on the AED when it arrives

6

Follow the voice instructions

7

Make sure the area is safe

1

Check the injured person's responsiveness

2

Call 999 immediately

3

Ask someone nearby to bring an AED

4

Note: If you're untrained, follow the emergency operator's and the device's instructions step by step



Turning On the AED

Turn on the device immediately and follow its instructions



**Press the
power button**



**Listen to the
instructions
attentively**



**Don't skip the
device's
instructions**



**Ask for calm
and for people
to stand clear**



**Continue until
the specialists
arrive**

Awareness Message: The device is designed to guide the user step by step — the most important action is turning it on and following its instructions



Preparing the Injured Person for AED Use

Prepare the chest area per the device's instructions



The injured person on a firm surface



Expose the chest area



Make sure the chest is dry



No pads over clothing



No water around the injured person



Secure the area before starting

Note: Don't delay requesting help or CPR while seeking perfection in preparation



Placing the AED Pads

The pads are placed on the chest so the device can analyze



Open the pad package

1

Place the pads according to the diagrams

2

Make sure they adhere well

3

Don't touch the adhesive surface

4

Use pediatric pads if available

5

Follow the device's instructions

6

Important Message: Follow the diagrams and guidance on the pads — they're designed to show the correct placement



The Device Analyzing the Case

Don't touch the injured person during analysis

Stay away from the injured person during analysis

1

Ask everyone not to touch them

2

Wait for the device's instructions

3

Don't press random buttons

4

Don't move the injured person during analysis


5

Important Note: Touching the injured person during analysis can affect the device's operation or hinder the procedure



When the Device Recommends a Shock

Make sure everyone is clear before delivering it

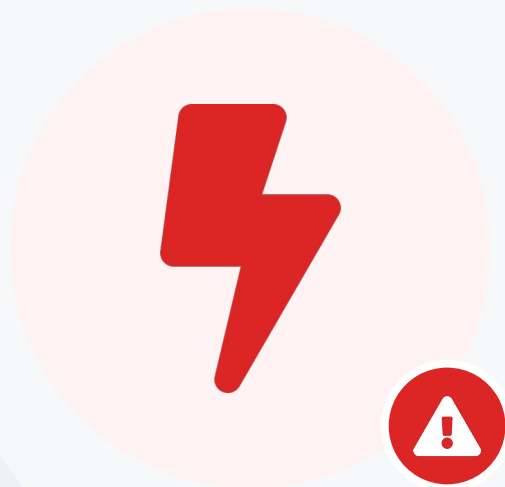

Say clearly: Stand clear

1

Visually confirm that no one is touching them

2

Follow the device's instructions precisely

3

Don't press the shock button without a prompt

4

After the shock, follow the CPR instructions

5

Watch for the device's next instructions

6

Important Message: The user doesn't decide to deliver the shock on their own — the decision is based on the device's analysis and instructions



If the Device Doesn't Recommend a Shock

No shock doesn't mean the case is minor



Don't deliver a shock without a recommendation



Follow the device's next instructions



Continue CPR if instructed



Monitor the injured person's responsiveness



Don't remove the pads without guidance

Note: You must keep requesting help and following the instructions until the specialists arrive



The Relationship Between AED and CPR

Two complementary skills in the initial response

AED



- Analyzing the heart rhythm •
- It guides when to stand clear and when to continue •
- It delivers a shock when needed •

CPR



- Temporarily supporting blood flow •
- Regular chest compressions •
- Don't stop without a reason •

Awareness Message: CPR and AED complement each other in the initial response to sudden cardiac arrest cases



Using an AED for Children

Special attention and following the device's instructions



**Pediatric pads
if available**



**Follow the
diagrams and
instructions**



**Follow the 999
operator's
instructions**



**Don't delay
requesting
help**



**Don't apply
what you
haven't been
trained in**

Important Note: Handling children requires approved hands-on training — request help immediately via 999



AEDs in Public Places

Know the device's location before an emergency occurs



Shopping malls



Airports



Schools and
universities



Sports
facilities



Hotels



Large
workplaces



Crowded public
facilities

Important Message: Knowing the AED's location before an emergency occurs shortens the response time when needed



The Role of Those Present at the Scene

Organizing roles reduces chaos



Someone calls
999



Someone
brings an AED



Someone
disperses the
onlookers



Someone
clears the way
for the
ambulance



Someone
follows the
device's
instructions



Someone
provides
information to
the specialists

Awareness Message: Organizing roles reduces chaos and helps with a quick response



Mistakes to Avoid

Actions that may delay the response or put those present at risk



**Delaying the
call to 999**



**Using the
device in a
dangerous
place**



**Touching the
injured person
during analysis**



**Touching them
while the shock
is delivered**



**Pressing the
shock button
without
guidance**



**Removing the
pads without a
reason**



**Crowding or
filming**



**Ignoring the
device's
instructions**

Important Message: Using an AED depends on calmness, following instructions, and not improvising randomly



Caring for the AED in Facilities

The device must be ready for use at any time



In a clear, easy location



Don't lock it away in a hard-to-reach place



A guidance sign for its location



Regular readiness checks



Validity of the pads and battery



Training staff in CPR and AED



Don't relocate it without documentation

Note: A device that isn't ready or whose location is unknown may not help during an emergency



After the Response Teams Arrive

Make room and provide the information



Move away if they ask



Tell them what happened



When the call was made



Whether CPR was performed



Whether an AED was used



Whether it recommended a shock



Follow their instructions

Important Message: After the specialists arrive, the role of those present becomes supporting organization and providing information

Do's and Don'ts When Using an AED

Correct Practices and Others to Avoid

Don't ✗

Don't delay requesting help

1

Don't use the device in a dangerous place

2

Don't touch the injured person during analysis

3

Don't press the shock button without guidance

4

Don't remove the pads without guidance

5

Don't crowd around the injured person

6

Don't film or post the case

7

Do ✓

Make sure the area is safe

1

Call 999 immediately

2

Ask for an AED to be brought

3

Turn on the device and follow its instructions

4

Stand clear during analysis and the shock

5

Continue until the specialists arrive

6

Learn CPR and AED from an accredited provider

7

Interactive Scenarios

Practical Application of What You've Learned

Scenario 2 An AED arrived but you're untrained

Correct Action: Turn on the device and follow the voice instructions, and stay on the line with 999

Scenario 1 A person in a mall collapsed and is unresponsive and not breathing normally

Correct Action: Secure the area, call 999, ask for an AED, and follow the device's and the emergency operator's instructions

Scenario 4 The device recommended a shock and people are nearby

Correct Action: Stop everyone, visually confirm they're clear, then follow the device's instructions precisely

Scenario 3 The device started analyzing and someone is touching the injured person

Correct Action: Move everyone away immediately, and don't allow them to be touched during analysis



Test Your Knowledge

Quick Questions to Assess Your Understanding

Q1: What is the first action when a person is unresponsive and not breathing normally?

Answer: Securing the area, calling 999, and asking for an AED to be brought if available



Q2: Is it acceptable to touch the injured person during AED analysis?

Answer: No, everyone must stay away from the injured person during the device's analysis



Q3: Does the user deliver the shock on their own?

Answer: No, the shock is delivered only if the device recommends it



Q4: Does an AED replace learning CPR?

Answer: No, the device is important but doesn't replace hands-on CPR training





Thank You

The AED is an effective tool when you know its location and follow its instructions precisely

Call 999 • Turn on the device • Stand clear during analysis and the shock • Follow the instructions precisely

"Knowing how to use an AED may make you the reason a life is saved"